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Lot number: 613100

Equipment delivery location:

mPTech Sp. z o.o.
Krakowska street 119
50-428 Wrocław
Poland

Designed in Poland
Manufactured in China

Manufacturer:

mPTech Sp. z o.o.
Nowogrodzka 31 street
00-511 Warsaw
Poland
Helpline: (+48 71) 71 77 400
E-mail: pomoc@mpotech.eu
Website: www.mpotech.eu

TABLE OF CONTENTS

1. Safety of use.....	10
2. Technical Specification	11
2.1 Technical data	11
2.2 Package content.....	11
3. Device construction.....	12
4. The use of the equipment.....	12
4.1 Charging the battery	12
4.2 Turning the device on, off and operating it.....	13
4.3 Replacing the suction nozzles	13
4.4 Use.....	13
5. Troubleshooting.....	14
6. Cleaning of the device	14
7. Correct disposal of used equipment and battery	15
8. The use of the manual	15
9. Warranty and service	16
10. Declaration of conformity with European Union Directives.....	16

1. Safety of use

Read these guidelines carefully. Not following them may be dangerous or illegal.

–INTENDED USE–

The device is intended for facial skin cleansing only and for private indoor domestic use only and must not be used for commercial purposes. Any other use is improper and dangerous. The manufacturer is not liable for injuries resulting from the use of the device in a manner inconsistent with its intended use. The device is not a medical equipment.

–CHILDREN SAFETY–

Keep the device, its accessories and packaging out of reach of children. Cables and plastic packaging may wrap around the child's neck or be swallowed, resulting in suffocation or choking.

Children must not play with this device, or clean it without an adult supervision.

–DO NOT RISK–

Do not use the device near flammable fabrics. Do not use the product with a damaged body or accessories. The product is not intended for use on animals. Do not use the device for purposes other than cleaning face.

–QUALIFIED SERVICE–

This product may only be repaired by qualified service personnel of the manufacturer or by an authorised service centre. Repairing the device by an unqualified or unauthorized service technician may result in damage to the device and void the warranty. The manufacturer is not responsible for problems caused by unauthorized modifications to the equipment made by the user.

–CHILDREN AND PHYSICALLY OR MENTALLY DISABLED PEOPLE–

This appliance may be used by children from 8 years of age and older and persons with limited physical, sensory or mental capabilities, or lack of experience and knowledge, if they are supervised or have been instructed in the safe use of the appliance and understand the hazards involved.

–WATER AND OTHER LIQUIDS–

Protect the device from water and other liquids. Avoid using the device in environments with high humidity like kitchens or bathrooms. Water particles may cause moisture to appear in the device, which adversely affects the electronics and may cause damage to the equipment.

–SURROUNDINGS–

The device should not be used in dusty conditions. While charging the battery make sure it cannot be damaged or endanger people or animals. Do not expose the device to sunlight and do not use it or leave it in close proximity to heat-generating devices (heaters, radiators etc.). Plastic parts may deform, which may affect the operation of the device or cause it to break down completely.

–BATTERY–

Avoid exposing the battery to very high/low temperatures (0°C – 45 °C / 32°F – 113°F). Extreme temperatures may influence capacity and battery life. Avoid contact with liquids and metal objects as it may cause partial or complete battery damage. The battery should be used according to its purpose. Do not destroy, impair or throw the battery into flames - it might be dangerous and may cause fire. Worn-out or impaired battery should be placed in a dedicated container. Excessive charging the battery can cause damage. Therefore do not charge longer than necessary. The battery in the device is not replaceable.

-CONTRAINDICATIONS-

Do not use the massager without a doctor's approval in the following cases: pregnancy, breastfeeding, skin cancer, presence of metal bone anastomoses in the areas where the device is used, fever, skin diseases, epilepsy, hypertension, heart failure, purulent dermatitis, anti-acne therapy with retinoids.

If you experience pain or discomfort while using the device, stop immediately.

2. Technical Specification**2.1 Technical data**

Facial cleansing device	
Model	BE-100124
Power	4 W
Suction power	65 kPa
Power supply	5 V = 500 mA
Battery	Li-ion 3,7 V = 500 mAh, 1,85 Wh
Additional	Four replaceable tips, micro USB cable, LCD

2.2 Package content

- Facial cleansing device
- Suction tips (4 pcs.)
- Filtering sponge (6 pcs.)
- Seals (2 pcs.)
- microUSB cable
- User Manual

If any of the above items is missing or damaged, please contact your dealer.


We advise to keep the container in case of any complaints.

You print a version of the user manual with the large font size or download it from the web address below:

<https://mptech.eu/qr/huslog>

Suction tips, filtering sponges and seals can be purchased from the manufacturer
Helpline: (071) 71 77 400 or pomoc@mptech.eu

3. Device construction

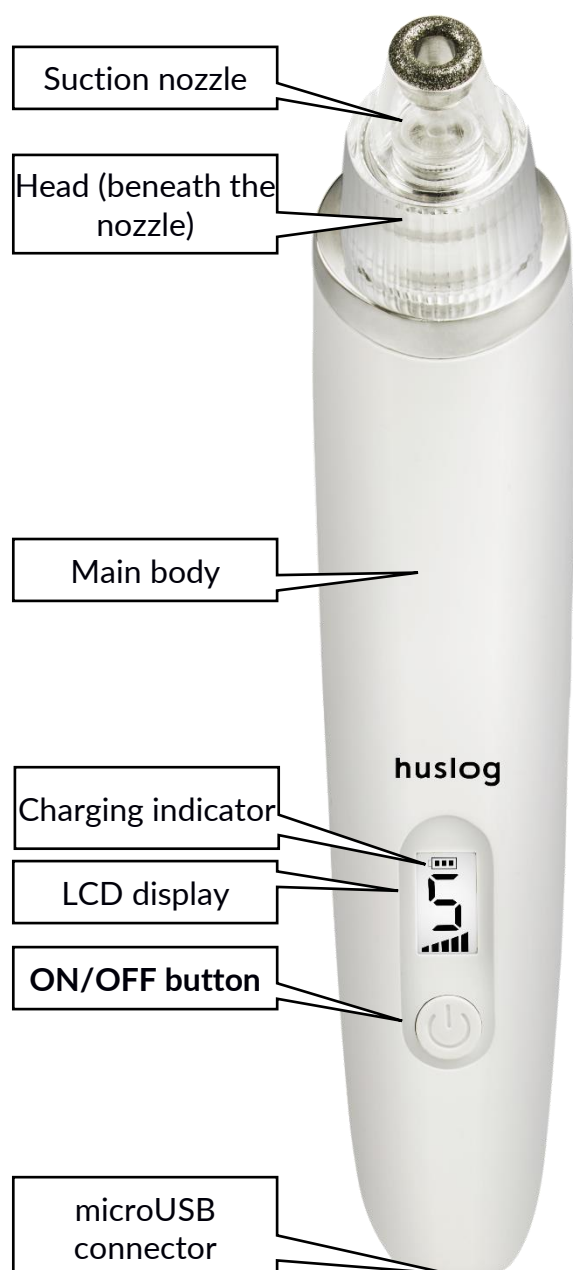The microcrystalline tip removes keratinised and dead skin and degreases the top layer of the skin.

The small round hole tip helps remove blackheads and whiteheads in hard to reach areas.

The tip with a large round hole gently massages, stimulates and tightens the skin.

The tip with an oval opening is designed for fine facial lines, enhancing skin texture and elasticity.

4. The use of the equipment

4.1 Charging the battery

 *The battery should be fully charged before the first use of product. Several charging and discharging cycles may be necessary to achieve full battery capacity.*

 *Avoid excessive battery discharge. This can damage the battery or shorten its service life. A low battery charge is indicated by a blank charge indicator symbol on the LCD display of the device.*

 *The device can heat up during charging. It is a regular occurrence and should not affect the stability or performance of the device. Do not leave the device to be charged unattended as this may result in fire, damage to the device or its surroundings etc.*

To charge your device, connect the microUSB plug of the supplied cable to the charging socket on your device. Connect the other end of the cable to a powered USB adapter or, for example, a socket on your computer, a powerbank, etc., whose output guarantees a power supply of 5 V = 500 A or higher. Once charging begins, the LCD display will illuminate for a few seconds and the charge indicator icon will flash according to the battery charge:

- Flashing first indicator segment - from 0% to approx. 33% of charging
- Flashing second indicator segment - from approx. 33 % to approx. 66% of charge.
- Flashing third indicator segment - from 66% of charge and above

When charging is complete, all segments of the charge indicator will be permanently displayed. Once charged disconnect USB from the device and from the power source.

4.2 Turning the device on, off and operating it

Press the ON/OFF button to turn the device on. The device will turn on at the lowest suction intensity. Continuously pressing the ON/OFF button will increase the suction intensity - from level 1 to level 5 (the current level is shown on the display of the device). When the maximum suction intensity is reached, another press of the ON/OFF button will switch the unit off. The device does not switch off automatically.

4.3 Replacing the suction nozzles

 *Only change tips when the machine is switched off.*

To remove the tip, grasp the body of the device with one hand and the corrugated part of the tip with the other. Pull the two pieces in opposite directions while rotating the tip around its axis.

The tip is installed by placing it on the head of the device and pressing it down. After installation, make sure that the tip is correctly positioned. Stable and secure assembly is ensured by two seals located on the head of the device.

After removing the tip it is possible to replace the filter sponge. Remove the device head plug and rotate the device with the head opening facing downwards - the filter sponge will fall out of the head. Replace the dirty filter sponge with a fresh one, block the head opening with a plug and attach the suction nozzle.

4.4 Use

 *Avoid using the device for too long and intensively in one place. It can cause skin irritation.*

 *Do not use the device on an irritated or damaged skin.*

 *The treated skin should be taut. Tighten the facial skin with your free hand by moving it down or to the side.*

REMOVAL OF OLD EPIDERMIS

Cleanse and dry your face. Place a warm towel on your face or perform a facial steam for a few minutes - this will open the pores of your face. Place the microcrystalline tip on the head of the device, turn it on and adjust the suction intensity. Then place the device on your

forehead and move from the centre outwards to the left and right side. Then place the device on your nose and move from the back of your nose towards your cheekbone to the left and right. Move to your chin and move the device from bottom to top. Finally, apply the device to the cheeks moving from the centre of the face outwards. Move from the bottom to the top of your face. After use, wash your face with soothing lotion and apply sunscreen.

SKIN CLEANING

Clean and dry your face. Place a warm towel on your face or perform a facial steam for a few minutes - this will open the pores of your face. Place the nozzle with the small or large round hole on the device head and switch on the device and adjust the suction intensity. Point the nozzle directly at the removal area. When removing blackheads, move the device slowly in one direction - do not make circular movements. For acne, apply the nozzle to the desired area and release immediately after suction. After use, rinse face with toner or use other pore-closing product.

REMOVAL OF FINE WRINKLES

Place the nozzle with the oval opening on the device head and switch on the device and adjust the suction intensity. Use the tip to gently massage the corners of the mouth, eyes and around the nose. Working with this handpiece helps improve blood circulation, skin elasticity and helps eliminate wrinkles

Immediately after using the device with any of the tips, the skin may become red and slightly irritated. Avoid chlorinated water, solariums, saunas and exposing your skin to the sun.

5. Troubleshooting

If you have problems with the unit, please refer to the following tips.

The device does not switch on	Check that the battery is charged. Check that the cables, plugs and sockets supplying power to the unit are not damaged.
After switching on, the suction does not work	Perhaps there has been a damage to the engine. Contact the manufacturer's service department.
If the problem persists, contact the manufacturer's service department.	

6. Cleaning of the device

To extend the service life of the device:

- **WARNING!** The device must not be submerged in water or washed under running water!
- For cleaning use only a dry cloth. Never use any agent with a high concentration of acids, alkalis or abrasives.
- Store the device and its accessories in a dry place away from the reach of children, sources of heat and moisture.
- Avoid very high temperatures, for it can shorten the life of electric and electronic components, melt the plastic parts and destroy the battery.

The packaging of the device can be recycled in accordance with local regulations.

Symbol for corrugated cardboard (packaging material).

PAP

The packaging elements should be disposed of in appropriate containers.

Marking of the container into which the stamping should be thrown - plastic/metal.

The marking of the container into which the box - paper should be discarded.

7. Correct disposal of used equipment and battery

The device is marked with a crossed-out garbage bin, in accordance with the European Directive 2012/19/EU on used electrical and electronic Equipment (Waste Electrical and Electronic equipment - WEEE). Products marked with this symbol should not be disposed of or dumped with other household waste after a period of

use. The user is obliged to get rid of used electrical and electronic equipment by delivering it to designated recycling point in which dangerous waste are being recycled. Collecting this type of waste in designated locations, and the actual process of their recovery contribute to the protection of natural resources. Proper disposal of used electrical and electronic equipment is beneficial to human health and environment. To obtain information on where and how to dispose of used electrical and electronic equipment in an environmentally friendly manner user should contact their local government office, the waste collection point, or point of sale, where the equipment was purchased.

According to the EU directive 2006/66/EC with changes contained in EU directive 2013/56/UE on the disposal of batteries, this product is marked with crossed-out garbage bin symbol. The symbol indicates that the batteries used in this product should not be disposed with regular household waste, but treated according to the

legal directives and local regulations. It is not allowed to dispose of accumulators with unsorted municipal waste. Users of battery and accumulators must use the available collection points of these items, which allows them to return, recycle and disposal. Within the EU, the collection and recycling of batteries and accumulators is a subject to separate procedures. To learn more about existing procedures of recycling of batteries and accumulators please contact your local office or an institution for waste disposal or landfill.

8. The use of the manual

Without prior written consent of the manufacturer, no part of this manual may be photocopied, reproduced, translated or transmitted in any form or by any means, electronic or mechanical, including photocopying, recording or storing in any systems of storing and

sharing information.

All brand names and product names are trademarks or registered trademarks of their respective companies.

9. Warranty and service

The product is covered by a 2-year warranty, starting from the date of sale of the product, except for accessories (USB cable, filter sponges, gaskets, suction tips) which are covered by a 6-month warranty. If you need to make use of a warranty claim, please contact the manufacturer's helpline or point of sale. The product returned for warranty claim should be complete and best in its original packaging.

10. Declaration of conformity with European Union Directives

mPTech sp. z o.o. hereby declares that the product complies with all relevant provisions of the so-called 'New Approach' of Directive of the European Union.